

What You Need

Item	Details
Beta Driver	Tobii.EyeTracker5.Offline.Installer_5.0.0.6035.exe (5.3 MB)
Download Link	<a href="#">TousLesDrivers</a>
Windows Version	Windows 10 (1709+) or Windows 11
Tobii Experience	Auto-installs from Microsoft Store (don't download separately)

What You Need

Item	Details
Beta Driver	Tobii.EyeTracker5.Offline.Installer_5.0.0.6035.exe (5.3 MB)
Download Link	<a href="#">TousLesDrivers</a>
Windows Version	Windows 10 (1709+) or Windows 11
Tobii Experience	Auto-installs from Microsoft Store (don't download separately)

Step 1 : Remove All existing Tobii Software : start fresh

You must remove all old tobii drivers and softwares first to be able to use tobii eye tracker 5 with windows eye control in presence of old tobii driver

- **Uninstall:** Go to **Settings > Apps & features**. Search for Tobii. Uninstall (tobii Driver, Experience App).
- It should look like this as whons in the image

Apps & features

App execution aliases

Search, sort, and filter by drive. If you would like to uninstall or move an app, select it from the list.

tobii

Sort by: Name Filter by: All drives

We couldn't find anything to show here. Double check your search criteria.

- **Check Hidden Files:** Press Win + R, type devmgmt.msc, and press Enter.
  - Expand **Software components**. If you see "Tobii," right-click and **Uninstall device**.
- **Restart:** Restart your computer.

Step 2: Install the Beta Driver (The Engine)

This specific driver adds the "HID" compatibility that Windows needs.


1. **Download:** Use the Beta Driver file: Tobii.EyeTracker5.Offline.Installer\_5.0.0.6035

Download Link [TousLesDrivers](#)

Download Link (Google Drive)

Please be careful while downloading , these are beta driver and sometimes they behave differently

Start installing

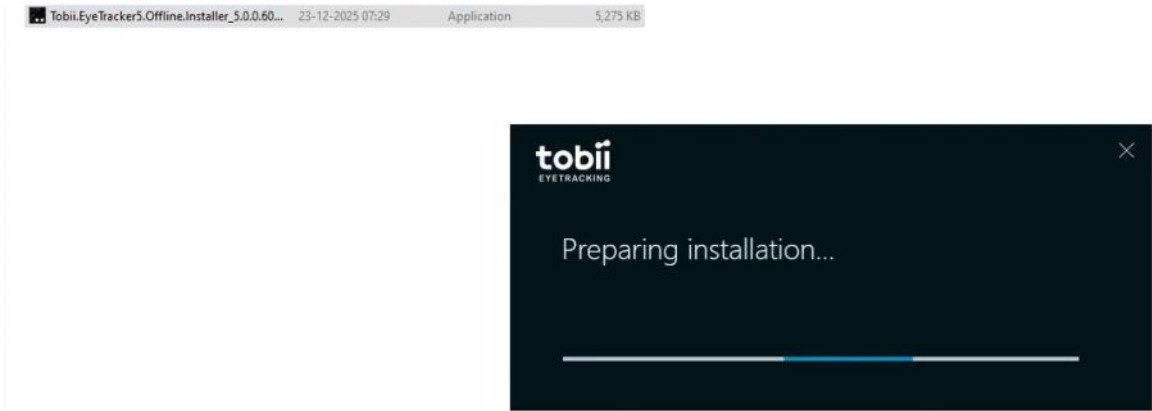

- 2 **Install:** Right-click the file and select **"Run as administrator"**. Follow the prompts.
- 3 **Restart:** Restart your computer immediately after installation.

## Apps & features

### App execution aliases

Search, sort, and filter by drive. If you would like to uninstall or move an app, select it from the list.

Sort by: Name ▾ Filter by: All drives ▾

1 app found




Tobii Device Drivers For Windows (EyeTracker5)

6.82 MB  
23-12-2025

### Step 3: Connect & Calibrate (The Steering)

Now we let Windows fetch the calibration app automatically.


1. **Plug In:** Connect your Tobii Eye Tracker 5 to a USB port.
2. **Wait:** Windows will detect the device and automatically download the **Tobii Experience App** from the Microsoft Store.
  - *Note: If nothing happens after 5 minutes, manually install "Tobii Experience" from the Microsoft Store.*
3. **Calibrate:** Open the **Tobii Experience** app. Create a profile and follow the dots to calibrate your eyes.

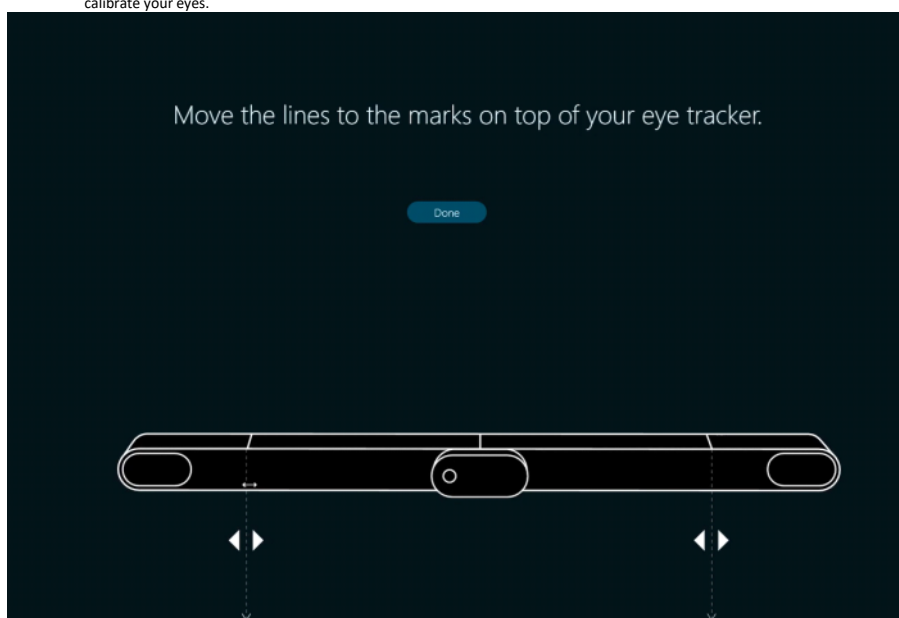


### Step 4: Enable Windows Eye Control (The Ignition)

1. Go to **Windows Settings > Accessibility (or Ease of Access)**.
2. Select **Eye Control**.

The "Launchpad" bar will appear at the top of your screen.

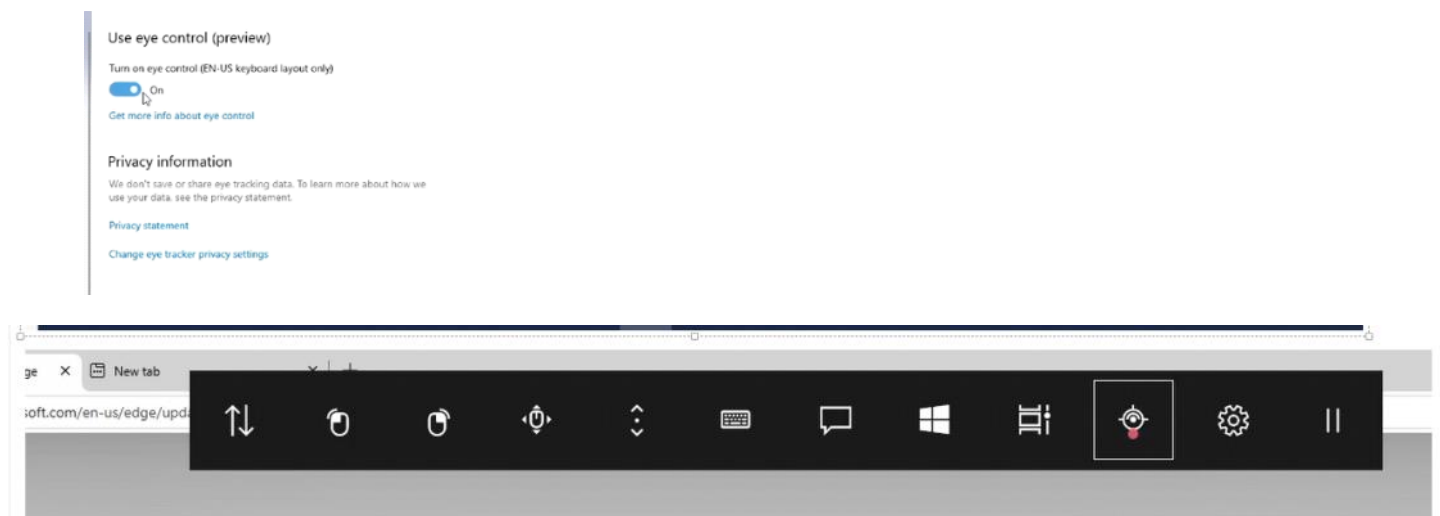

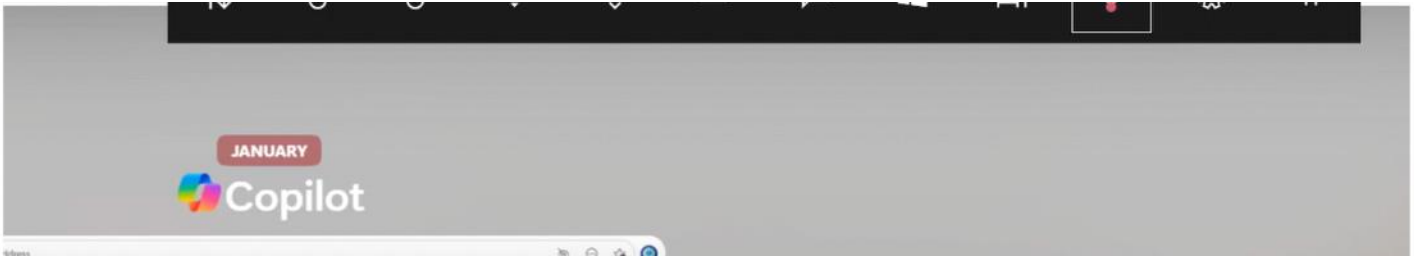

### Final Verification Checklist

Your setup is complete and working if:

- Beta driver 5.0.0.6035 shows in Apps & features
- Tobii Experience app installed and calibrated
- Eye tracker LEDs are visible when plugged in
- Eye Control toggle stays ON
- White Launchpad bar visible at top of screen
- Can see icons on Launchpad
- Looking at icons makes them highlight
- Dwelling activates the icon
- Can click, type, and control mouse with eyes

### What You've Installed (Summary)

After successful setup, you have:

Component	Version	Purpose
Beta Driver	5.0.0.6035	HID interface for Windows Eye Control
Tobii Experience	Latest from MS Store	Calibration and eye tracking
Windows Eye Control	Built into Windows	User interface (Launchpad)