

# Tobii Eye Tracker 5 with Windows Eye Control (Beta)

**Experimental setup** using a discontinued beta driver to enable Windows built-in Eye Control.

## About Tobii Eye Tracker 5:

- Gaming device — NOT designed for accessibility
- Does NOT support Windows Eye Control natively
- Requires **beta driver installation** to use Windows Eye Control
- Beta driver is discontinued — may have stability issues

## What You Need

### HARDWARE

[Tobii Eye Tracker 5 — Amazon India](https://www.amazon.in/s?k=tobii+eye+tracker+5) (<https://www.amazon.in/s?k=tobii+eye+tracker+5>)

### SOFTWARE DOWNLOADS

[Beta Driver v5.0.0.6035 \(Community Archive\)](https://www.touslesdrivers.com/index.php?v_page=23&v_code=69377) ([https://www.touslesdrivers.com/index.php?v\\_page=23&v\\_code=69377](https://www.touslesdrivers.com/index.php?v_page=23&v_code=69377)).

Warning: This uses a discontinued beta driver with no official Tobii support.

## System Requirements

- **Windows 10 or 11** (64-bit)
- **RAM:** 8 GB or more



## Quick Steps Overview

- 1 Remove old Tobii files
- 2 Install Beta Driver v5.0.0.6035
- 3 Verify in Apps & Features
- 4 Connect Tobii Eye Tracker 5
- 5 Display Setup & Calibration in Tobii app
- 6 Enable Windows Eye Control

1

### Remove Old Tobii Files (Clean Start)

Before installing the beta driver, remove old or conflicting Tobii files to ensure a clean installation:

#### 1. Check Apps & Features → search "Tobii"

- if no Tobii software or driver found, proceed to Step 2
- If found, uninstall **ALL** entries:
  - Tobii Experience Software
  - Tobii Eye Tracker 5 Drivers
  - Tobii Service
  - Any other Tobii entries

#### 2. Restart Computer

- Restart after uninstalling all Tobii software
- This ensures all services are fully stopped

#### Apps & features

##### App execution aliases

Search, sort, and filter by drive. If you would like to uninstall or move an app, select it from the list.

Sort by: Name ▾

Filter by: All drives ▾

We couldn't find anything to show here. Double check your search criteria.


*After cleanup: "We couldn't find anything to show here"*

## Download & Install Beta Driver

- Download the beta driver from: [TousLesDrivers](https://www.touslesdrivers.com/index.php?v_page=23&v_code=69377) ([https://www.touslesdrivers.com/index.php?v\\_page=23&v\\_code=69377](https://www.touslesdrivers.com/index.php?v_page=23&v_code=69377)) (warning: This is NOT an official Tobii download)
- File name: `Tobii.EyeTracker5.Offline.Installer_5.0.0.6035.exe`
- File size: 5.3 MB
- Right-click the installer and select "**Run as Administrator**"

### Start installing

Tobii.EyeTracker5.Offline.Installer\_5.0.0.60... 23-12-2025 07:29 Application 5,275 KB

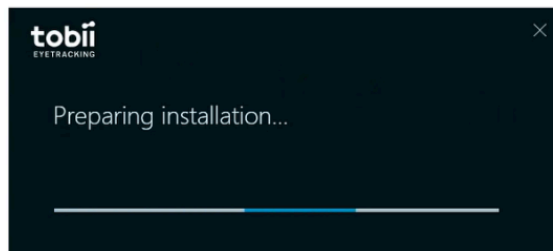


*Beta driver installation in progress*

3

## Verify Installation in Apps & Features

After installation completes, verify the driver was installed:

- Open **Settings** → **Apps & Features**
- Search for "**tobii**" again
- You should now see: "**Tobii Device Drivers For Windows (EyeTracker5)**"
- Size: ~6.82 MB

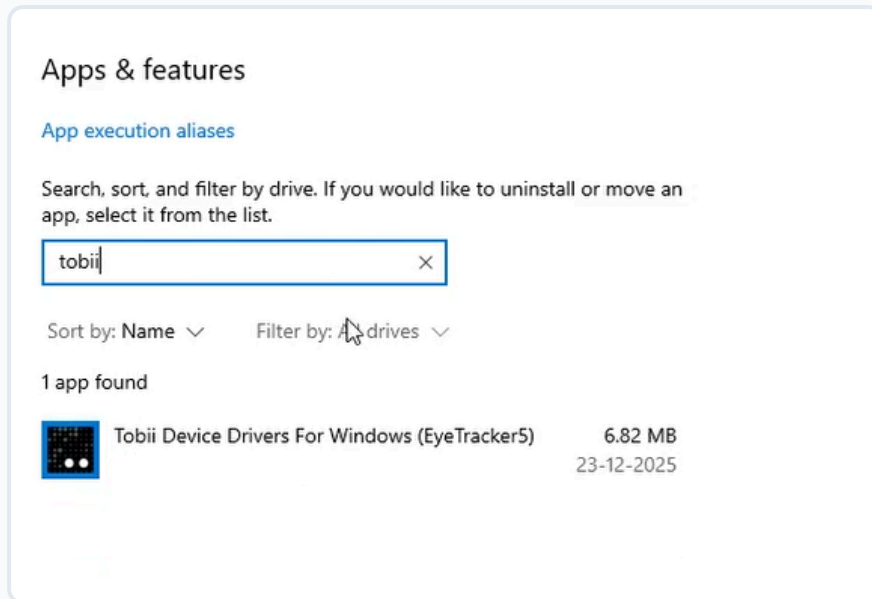

*Tobii Device Drivers successfully installed*

4


## Connect Tobii Eye Tracker 5

- Attach the Eye Tracker 5 to the bottom center of your monitor using the magnetic mount or double-sided tape
- Plug the Eye Tracker 5 into a USB port on your computer

## Display Setup & Calibration in Tobii app

### 1. Launch Tobii Experience & Set Up Display:


- Launch **Tobii Experience** from Start Menu (auto-installed from Microsoft Store)
- Click "**Set up display**" button when prompted
- Align the on-screen lines to the physical marks on top of your eye tracker

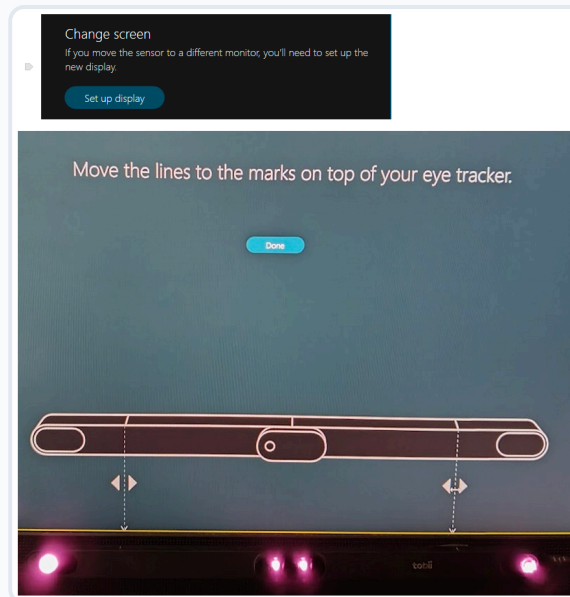

**Display Setup:** Click "Set up display" → Align the dashed lines with the physical marks on your Tobii Eye Tracker 5

### 2. Run Calibration:

- Position yourself within the tracking zone (45-95 cm from screen)
- Click "**Improve calibration**" in Tobii Experience
- Adjust the screen angle or slightly elevate PALS bed until the two white dots are visible in the center

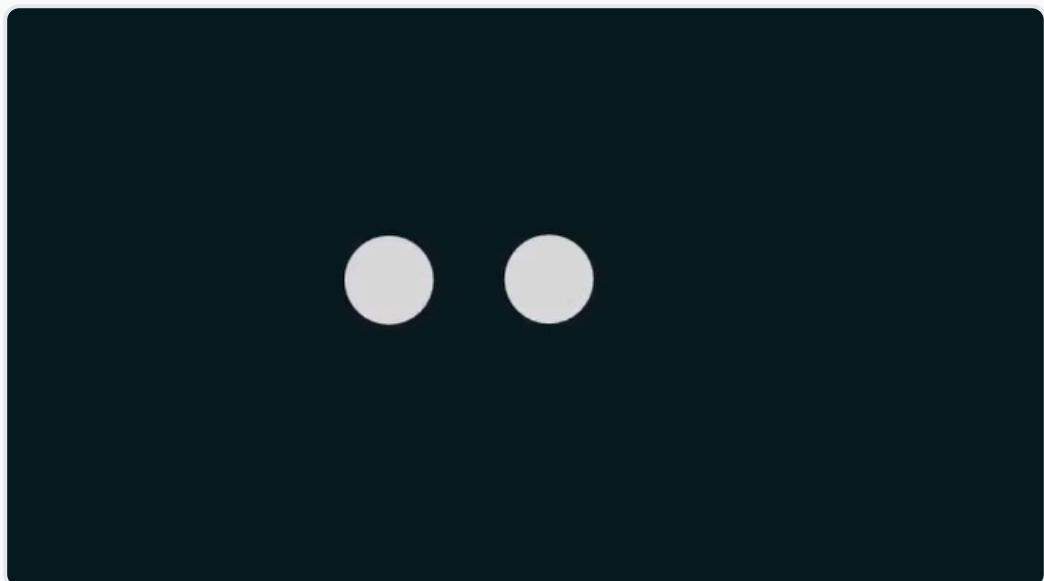

- Follow dots on screen with **EYES ONLY** — do not move your head
- Look at each dot until it disappears




**Calibration Screen:** Look at each blue dot until it disappears — keep your head still

## Enable Windows Eye Control

### Navigate to Windows Eye Control Settings:

- Open **Settings** → **Ease of Access** (or **Accessibility** in Windows 11)
- Select **Eye Control** from the left menu (under "Vision" section)
- Toggle "**Turn on eye control (EN-US keyboard layout only)**" to **ON**
- **RESTART your computer** if windows eye control doesn't work
- After restart, Eye Control should recognize the tracker
- The Eye Control launchpad (toolbar) will appear at the top of screen

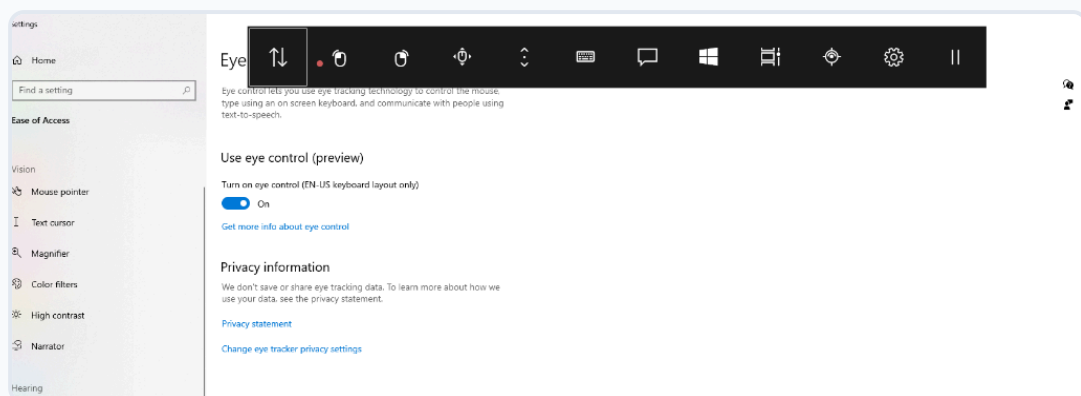

*Windows Settings → Ease of Access → Eye Control with toggle ON*

**Important:** You may need restart your computer after enabling Eye Control. The feature may not work properly without a restart sometimes.

### Known Issues with Beta Driver:

- May crash on startup with newer Windows builds
- Compatibility issues with some Windows updates
- No official support from Tobii — driver is discontinued

**If you experience problems:** Switch to Tab 1 (OptiKey) for a stable experience.

**Setup Complete!** You now have: Beta Driver v5.0.0.6035 and Windows Eye Control enabled. Use the Eye Control launchpad for keyboard, mouse, and navigation.