

# Tobii Eye Tracker 5 with OptiKey Setup

## About Tobii Eye Tracker 5:

- Gaming device — NOT designed for accessibility
- Does NOT support Windows Eye Control natively (Requires beta driver installation to use Windows Eye Control)
- compatible with **OptiKey**

## What You Need

### HARDWARE

[Tobii Eye Tracker 5 — Amazon India](https://www.amazon.in/s?k=tobii+eye+tracker+5) (<https://www.amazon.in/s?k=tobii+eye+tracker+5>).

### SOFTWARE DOWNLOADS

[Tobii Experience Driver v4.183.0.30025](https://gaming.tobii.com/getstarted/) (<https://gaming.tobii.com/getstarted/>).

[OptiKey Pro Classic v3.2.5](https://optikey.org/help/versions) (<https://optikey.org/help/versions>).

## System Requirements

- **Windows 10 or 11** (64-bit)
- **RAM:** 8 GB or more



## Quick Steps Overview

- 1 Install Tobii Experience Driver v4.183.0.30025
- 2 Connect Eye Tracker 5 hardware
- 3 Setup Tobii Experience app
- 4 Verify installation in Apps & Features
- 5 Install OptiKey Pro Classic v3.2.5
- 6 Use OptiKey for eye tracking

**Important:** OptiKey v4 does NOT support Eye Tracker 5. You must use **OptiKey Classic v3.2.5**.


1

## Download & Install Tobii Driver

This driver allows your computer to communicate with the eye tracker.

### 1. Download the Driver:


- Go to: [gaming.tobii.com/getstarted](https://gaming.tobii.com/getstarted/) (<https://gaming.tobii.com/getstarted/>)
- Under "1. Select the hardware", choose **"Tobii" (Eye Tracking)**
- Under "2. Select Tobii device", choose **"Tobii Eye Tracker 5"**
- Click **"DOWNLOAD DRIVER"**

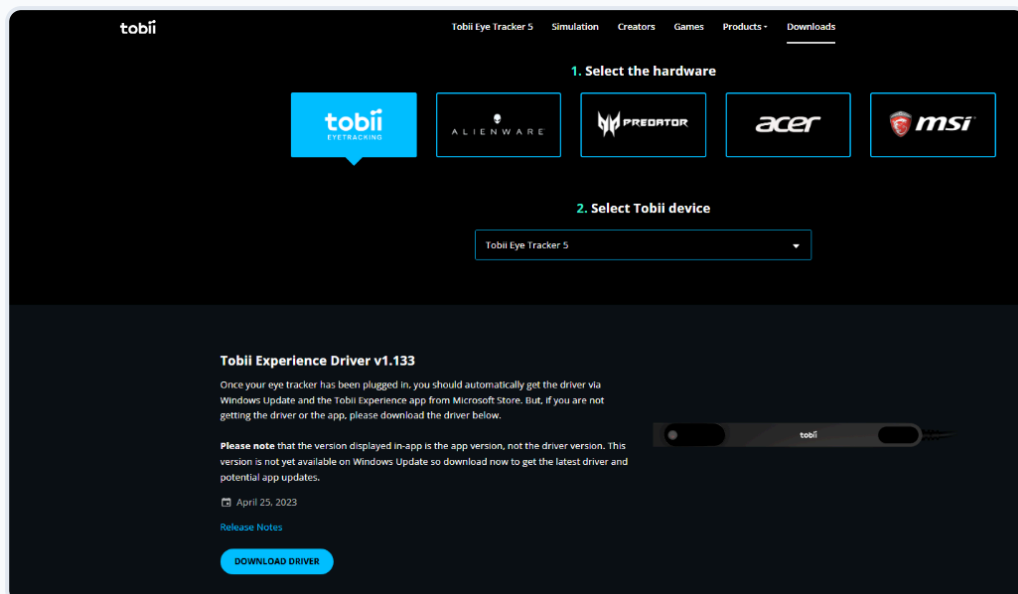

Select "Tobii Eye Tracker 5" and click "DOWNLOAD DRIVER"

### 2. Install the Driver:

- Find the downloaded file: `Tobii.EyeTracker5.Offline.Installer_4.183.0.30025`
- **Right-click** the file Select **"Run as administrator"**
- Follow the installation prompts, click "Next" until finished
- The installer will also install **Tobii Experience Software**

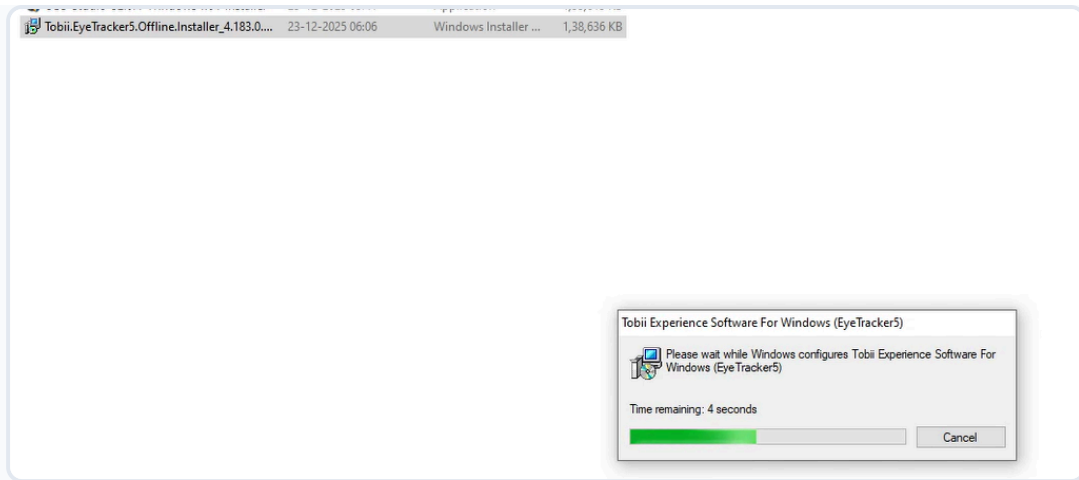

*Tobii Experience Software installing — wait for completion*

**After installation:** Check Apps & Features → search "Tobii" to verify the driver is installed

2

## Connect Tobii Eye Tracker 5


- Attach the Eye Tracker 5 to the bottom center of your monitor using the magnetic mount or double-sided tape
- Plug the Eye Tracker 5 into a USB port on your computer

## Setup Tobii Experience App

This app calibrates your eyes for tracking:

### 1. Open the App:


- Click the **Start Button**
- Type: **Tobii Experience**
- Open the app
- If not found: Open **Microsoft Store**, search "Tobii Experience", and install it

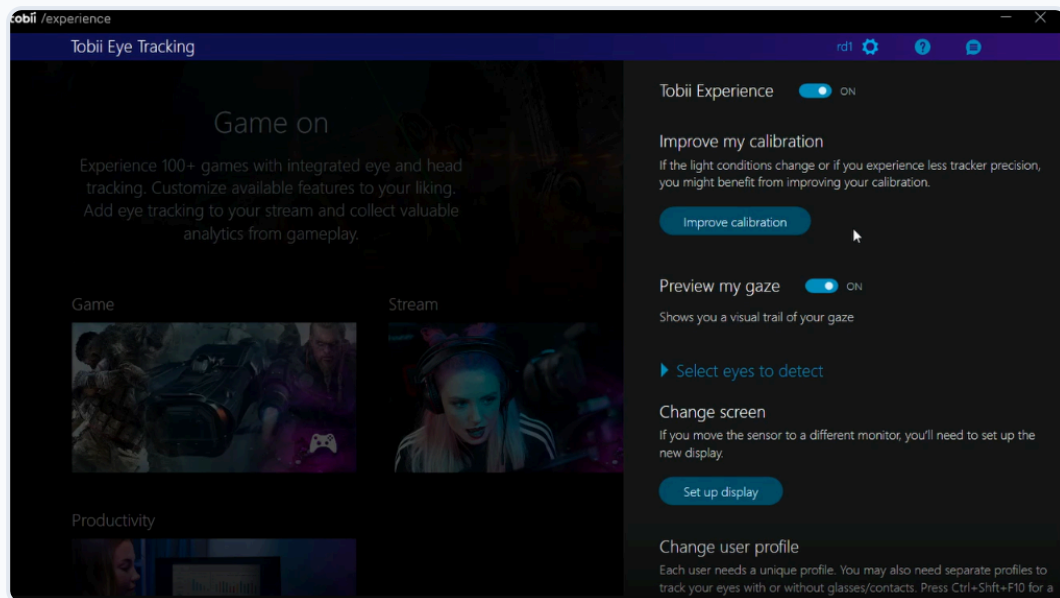

***Tobii Experience App:** Main interface with calibration and display setup options*

### 2. Display Setup:


- Click "**Set up display**" button
- Align the on-screen lines to the physical marks on top of your eye tracker
- Click "**Done**" when aligned

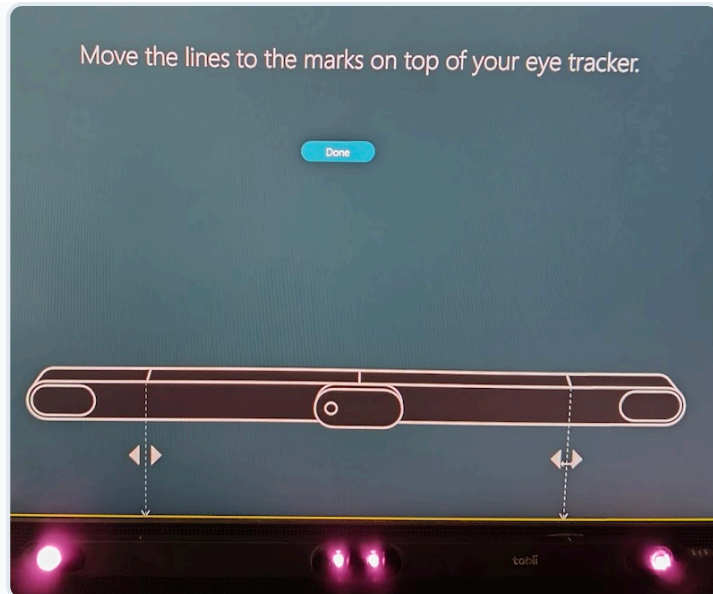

**Display Setup:** Move the dashed lines to align with the marks on your Tobii Eye Tracker 5

### 3. Calibrate:

- Click "**Improve calibration**" button
- Follow the dots on screen with **EYES ONLY** — keep your head still
- Look at each dot until it disappears
- Save your calibration profile when complete

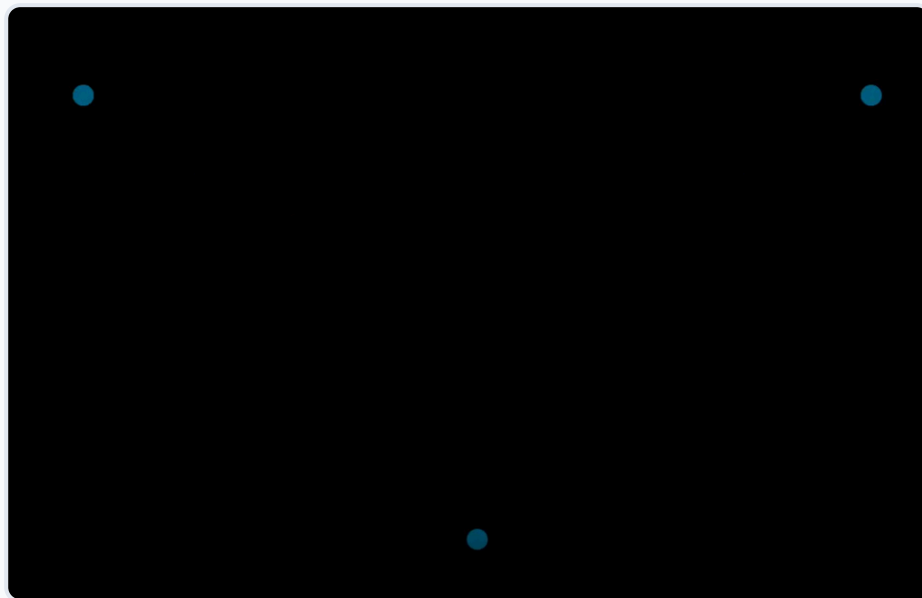


**Calibration:** Look at each blue dot until it disappears — keep your head still

**About Tobii Experience:** Use this app to improve calibration, recalibrate whenever needed, change profiles, and re-run Display Setup if you move the tracker.

## Verify Installation

After installation, check that all required software is installed:

- Open **Settings** → **Apps & Features**
- Search for **Tobii**
- You should see **3 apps**:
  - **Tobii Experience** (the app from Microsoft Store)
  - **Tobii Experience Software For Windows**
  - **Tobii Experience Software For Windows (EyeTracker5)** — the driver

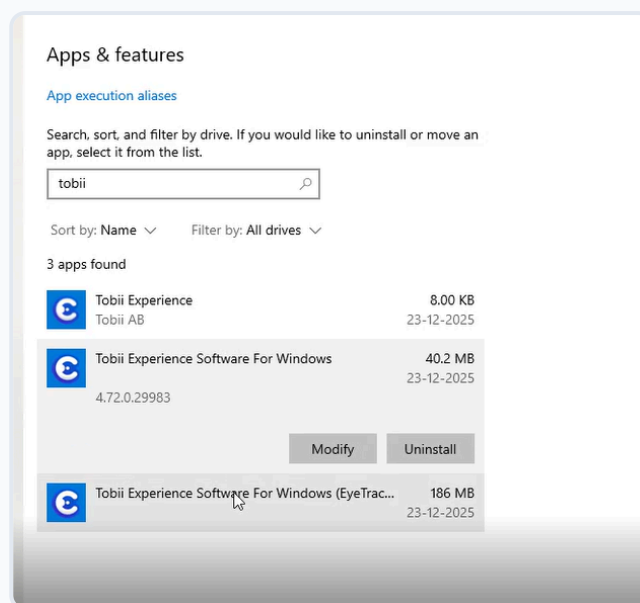

*Apps & Features should show these 3 Tobii components after successful installation*


## Install OptiKey

OptiKey is the actual software used for communication.

Do NOT download Version 4. It will NOT work with Tobii Eye Tracker 5.  
You **MUST** use **Version 3 (OptiKey Classic)**.

### 1. Download:


- Go to: [optikey.org/help/versions](https://optikey.org/help/versions) (<https://optikey.org/help/versions>)
- Scroll down to find the "OptiKey Classic (V3)" section
- Click: "Download OptiKey Pro Classic (v3.2.5)"

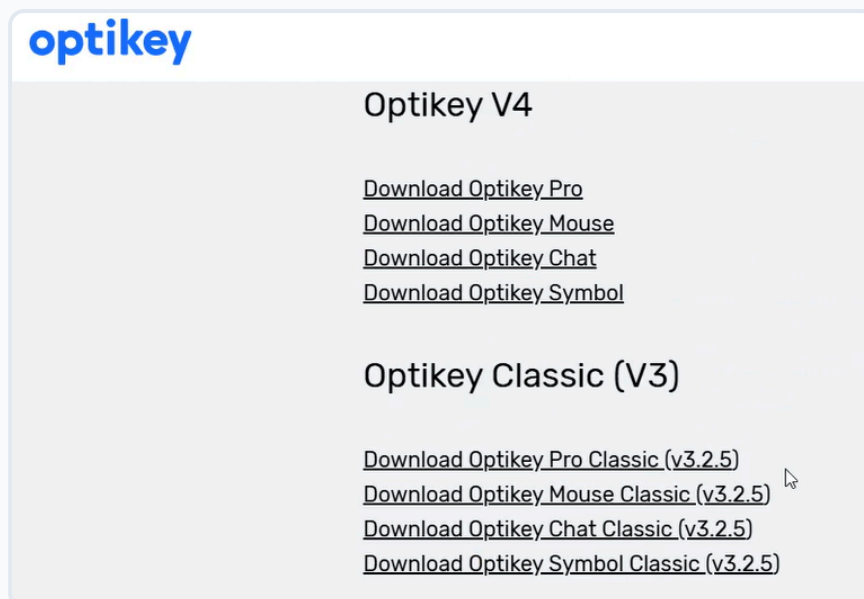

*OptiKey versions page — download from "OptiKey Classic (V3)" section, NOT V4!*

### 2. Run the Installer:


- Find the downloaded file: OptiKeySetup-3.2.5-Pro.exe
- **Right-click** → "Run as administrator"
- The Setup Wizard will start

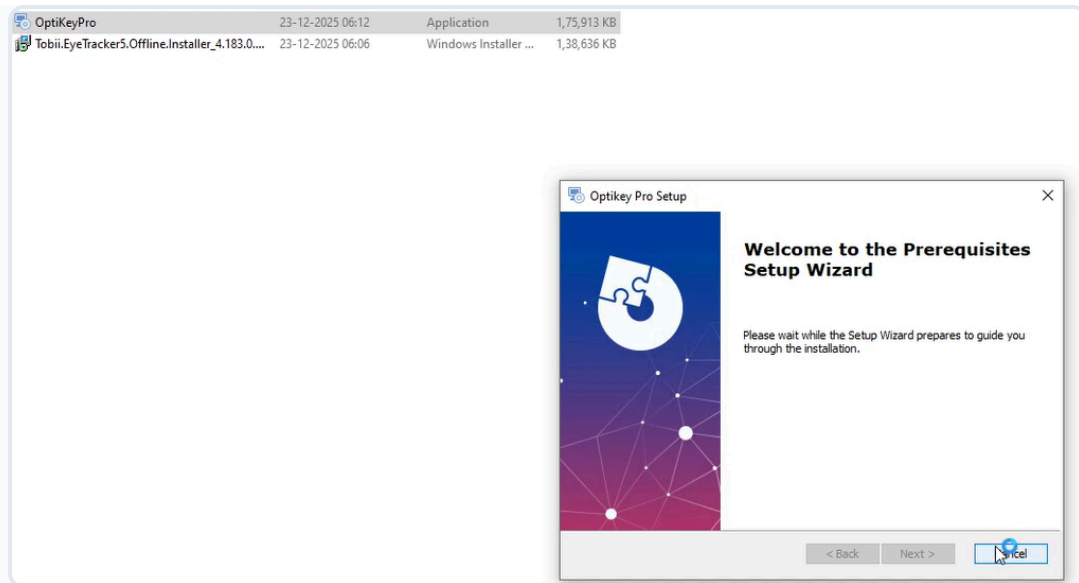

*OptiKey Pro Setup Wizard — click "Next" to continue*

### 3. Select Eye Tracker:

- When prompted, select "**Tobii Eye Tracker 5**" from the dropdown
- This configures OptiKey to work with your device

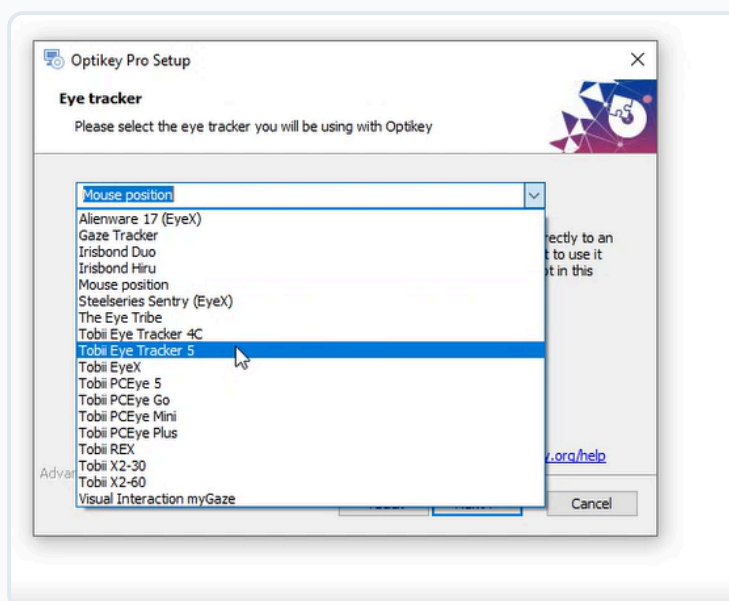


*Select "Tobii Eye Tracker 5" from the eye tracker dropdown*

### 4. Complete Installation:

- Click "**Install**" to begin the installation
- Wait for the installation to complete



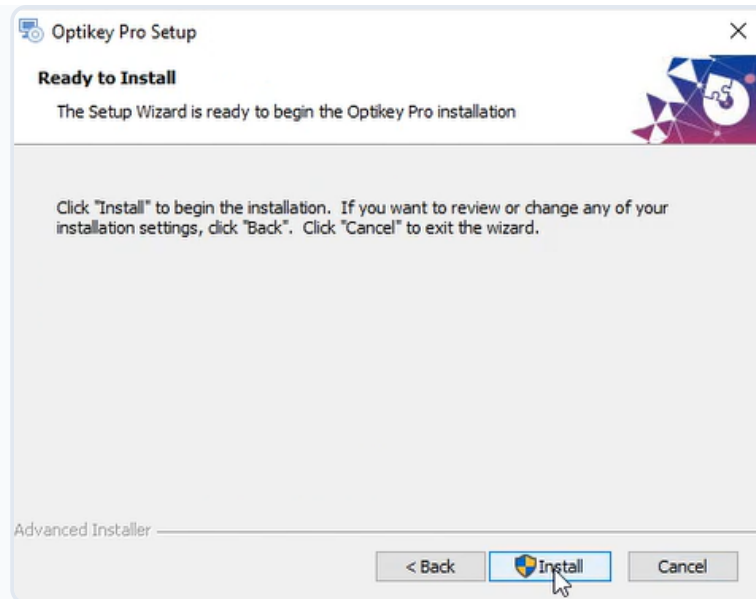

*Click "Install" to complete the OptiKey Pro installation*

## Use OptiKey for Eye Tracker 5

### 1. Open OptiKey:


- Launch the program from the **Start Menu**
- The OptiKey keyboard will open on your screen

### 2. Test & Verify:

- The keyboard should show **"Pointing: Tobii Eye Tracker 5"** in the info panel
- Look at keys to type — your gaze controls the cursor
- Try typing a message to test the eye tracking



*OptiKey Pro v3.2.5.0 showing "Pointing: Tobii Eye Tracker 5" — Setup confirmed!*

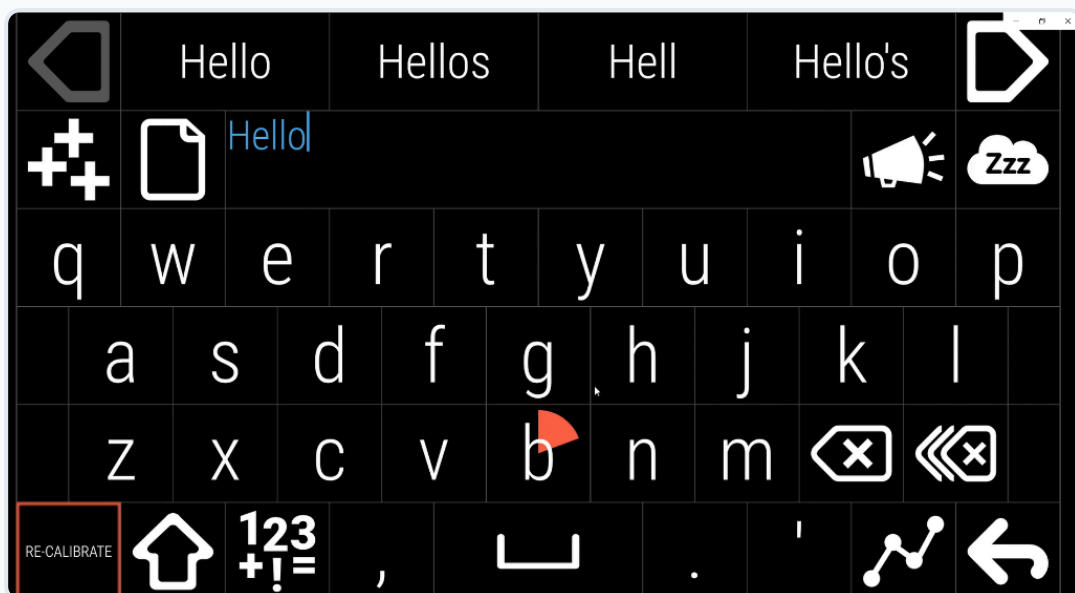

*OptiKey keyboard in action — type by looking at keys*

### Tips & Settings:

- **Recalibrate:** Click "RE-CALIBRATE" button on OptiKey keyboard (bottom-left)
- **Dwell Time:** Adjust in Settings → Pointing & Selecting
- **Language:** Hindi keyboard available! Go to Settings → Words



*OptiKey Pro Classic v3.2.5 working with Tobii Eye Tracker 5*

**Setup Complete!** You now have: Tobii Experience Driver v4.183.0.30025, Tobii Experience app, and OptiKey Pro Classic v3.2.5. Full eye tracking control of keyboard, mouse, and speech output is now available through OptiKey.