

Tobii Eye Tracker 4C with Windows Eye Control Setup

Good option for Windows Eye Control — used by many PALS.

Discontinued: Tobii no longer sells the Eye Tracker 4C. You will need to purchase a used unit from eBay or other second-hand marketplaces.

About Tobii Eye Tracker 4C:

- **Native Windows Eye Control support** — no workarounds needed
- Discontinued — only available second-hand now

What You Need

HARDWARE

[Tobii Eye Tracker 4C — Available on eBay \(used\)](https://www.ebay.com/sch/i.html?_nkw=tobii+eye+tracker+4c) (https://www.ebay.com/sch/i.html?_nkw=tobii+eye+tracker+4c)

eBay — ₹15,000-20,000 (used)

SOFTWARE DOWNLOADS

[Tobii Core Software \(for 4C\)](https://gaming.tobii.com/getstarted/?bundle=tobii-eye-tracker-4c) (https://gaming.tobii.com/getstarted/?bundle=tobii-eye-tracker-4c)

System Requirements

- **Windows 10** (64-bit) — Core Software is NOT compatible with Windows 11
- **RAM:** 8 GB or more



Quick Steps Overview


- 1 Download from gaming.tobii.com/getstarted (https://gaming.tobii.com/getstarted/) → Select 4C & Core Software
- 2 Install software
- 3 Verify in Apps & Features → "Tobii Eye Tracking" installed
- 4 Connect Tobii Eye Tracker 4C
- 5 Display Setup & Calibration in Tobii app
- 6 Enable Windows Eye Control in Settings

Step-by-Step Installation Guide

1

Download Tobii Core Software

- Go to: gaming.tobii.com/getstarted/ (<https://gaming.tobii.com/getstarted/>).
- Select "Tobii Eye Tracker 4C"
- Select "Tobii Eye Tracking Core Software"
- Download **Tobii Eye Tracking Core Software v2.16.8**

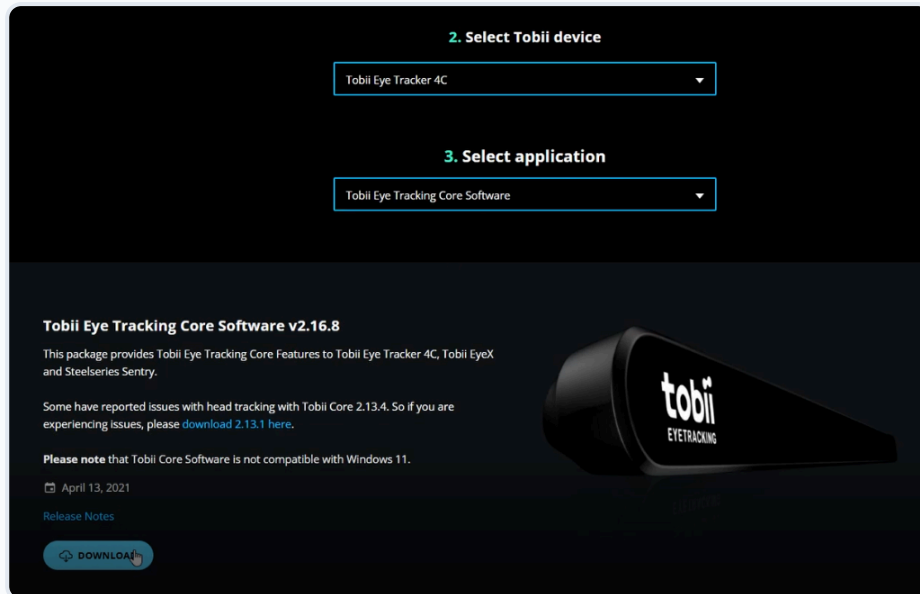


Download Page: Select "Tobii Eye Tracker 4C" and "Tobii Eye Tracking Core Software", then click Download

2

Install & Restart

- Run the downloaded installer
- Complete installation
- Restart your computer (optional)

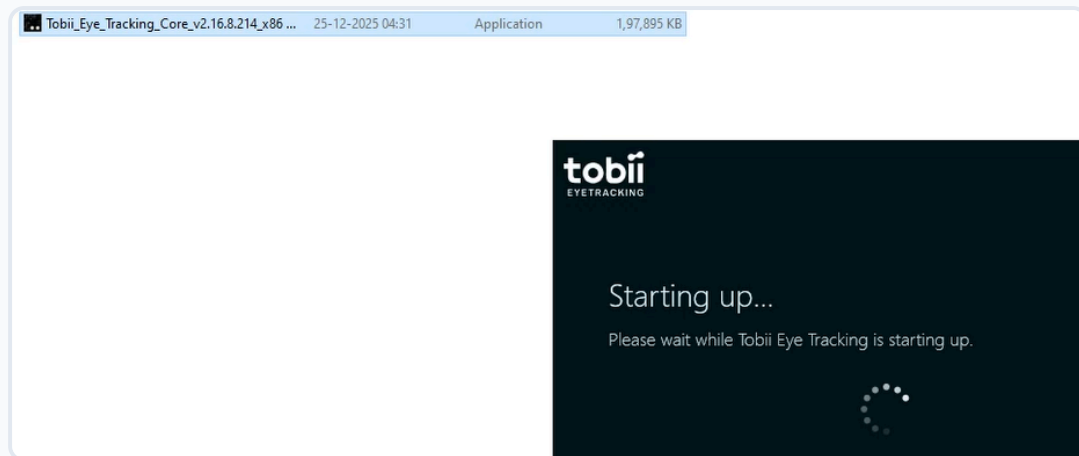


Installation: Run the downloaded installer and wait for Tobii Eye Tracking to start up

3

Verify Installation

- Check settings → Apps & Features → search "Tobii" → should show v2.16.8

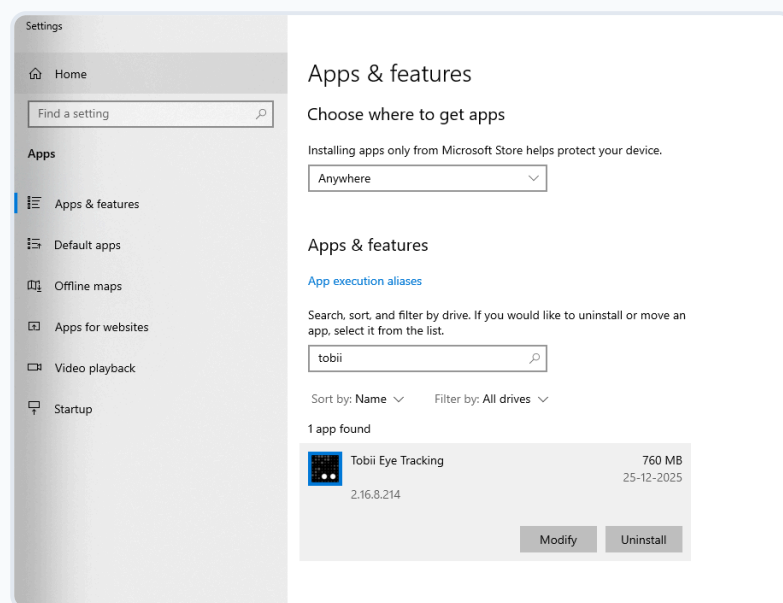

Verify Installation: Apps & Features should show "Tobii Eye Tracking" v2.16.8.214 (760 MB)

Connect Tobii Eye Tracker 4C


- Attach the Tobii 4C to the bottom center of your monitor using the magnetic strip or double-sided tape
- Plug into USB 3.0 port
- Wait for Windows to detect the device

Verify Tobii Service is Running:

- Press `Win + R` , type `services.msc`
- Find "Tobii Service" → should show "Running"

Display Setup & Calibration


- Launch Tobii Eye Tracking → Click "**Set up display**"
- Align the on-screen lines to the physical marks on top of your eye tracker

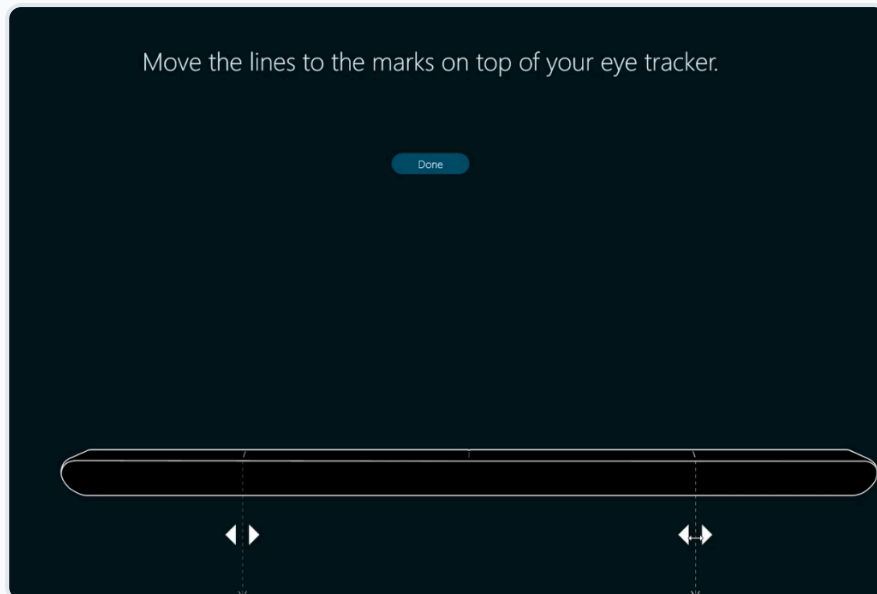

Display Setup: Move the dashed lines to align with the marks on your Tobii 4C

Run Calibration:


- Click "**Improve calibration**" in Tobii Eye Tracking app
- Look at each blue dot until it explodes — keep your head still
- Save your calibration profile after success

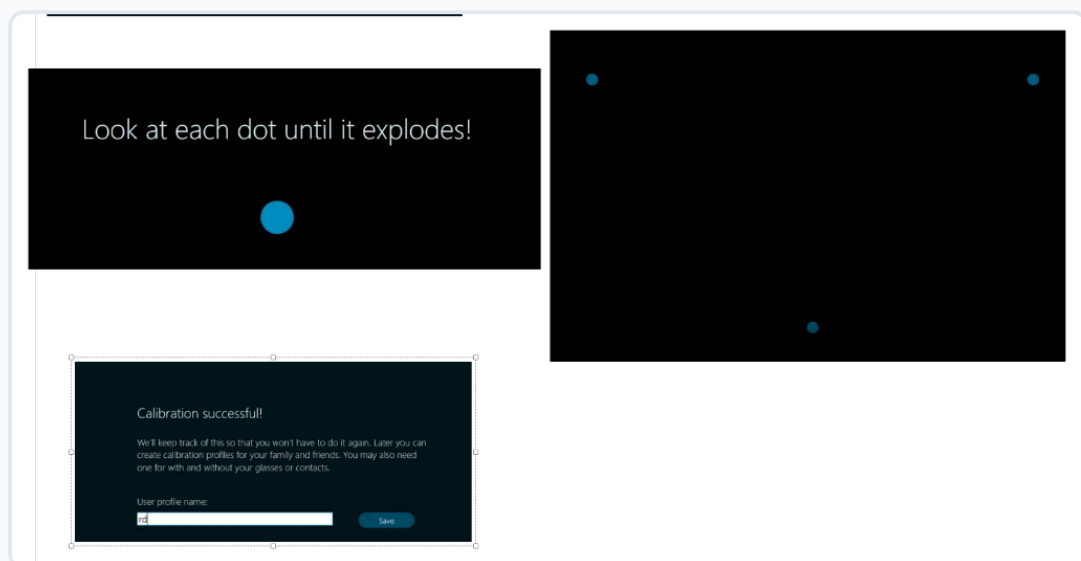

Calibration: Look at each blue dot until it explodes, then save your profile

Enable Windows Eye Control

- Go to: Settings → Ease of Access → Eye Control
- Toggle "Use eye control" to ON
- Eye Control toolbar will appear on screen

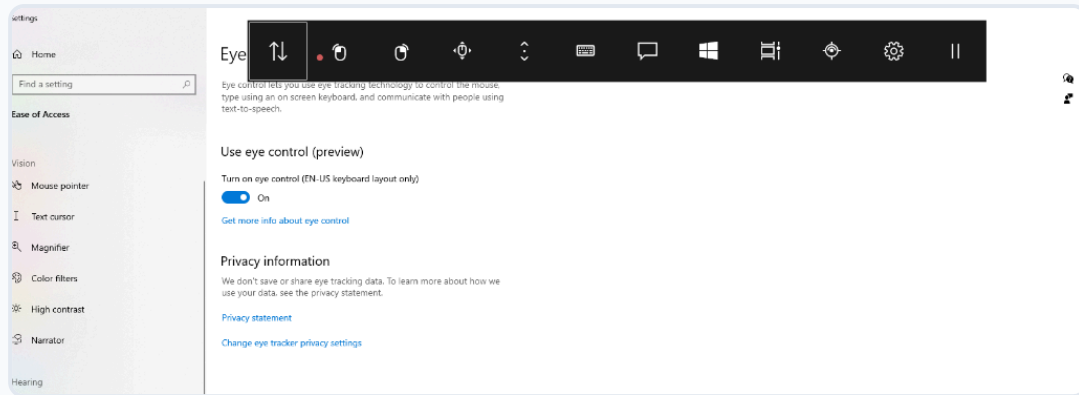

Enable Windows Eye Control from Ease of Access settings



Download Links

→ [Official Tobii Download Page](https://gaming.tobii.com/getstarted/) (https://gaming.tobii.com/getstarted/)

→ [Tobii Help Center](https://help.tobii.com/) (https://help.tobii.com/)